

Disposable lifting slings

SystemRoMedic™

Personal transfer aids for short-term use provide optimal hygiene and safety

Handicare's SystemRoMedic™ range offers an increasing selection of personal transfer aids for short-term safety use and optimal hygiene, among others a number of lifting sling models for seated or horizontal lifting as well as for sit-to-stand and standing training. Most of these lifting slings are made from a sturdy and durable non-woven (polypropylene) material but we also use a soft and durable cotton/polyester blend. Instead of laundering, the disposable lifting slings are discarded after use, when they become soiled or when the users no longer need them. This prevents the spread of multi-resistant and other types of infectious bacteria.




StretcherSling[™]

SystemRoMedic[™]

For horizontal lifting with StretcherBar

StretcherSling is a lifting sling that is used in combination with mobile lifts or stationary lift systems and Handicare's SystemRoMedic[™] StretcherBar, a specially designed sling bar for horizontal lifting and transferring of users in a lying position.

StretcherSling has ten suspension points; five on either side of the user, providing evenly distributed support under the entire body.

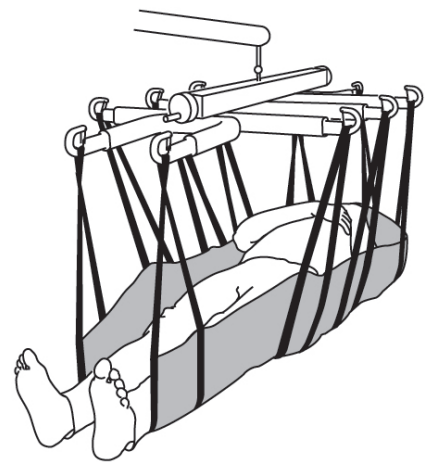

Item No.	Product	Size	Unit(s)
46592007	StretcherSling, disposable, non-woven	one size	5 pcs
	SWL 300 kg / 660 lbs		



SafeHandlingSheet[™]

SystemRoMedic[™]

For horizontal lifting and positioning with a two-point sling bar

Handicare's disposable SafeHandlingSheet is used for horizontal lifting or for positioning of a user in bed without subjecting the caregiver to unnecessary strain. SafeHandlingSheet is used in combination with mobile lifts or stationary lift systems and a two-point sling bar.

The disposable SafeHandlingSheet is made from a soft and breathable cotton/polyester blend. The two largest sizes are designed without straps running across the sheet which can increase pressure.

- For supine transfers, turning and repositioning
- Supports prevention of pressure sores
- Skin safe, breathable cotton blend material provides improved comfort for the user
- Infection control (use with patients on contact, droplet and airborne precautions)
- Compatible with lifts and sling bars from other manufacturers
- Bariatric options are available



Item No.	Product	Size	Unit(s)
47512006	SafeHandlingSheet, disposable, cotton	M	5 pcs
47512007	SafeHandlingSheet, disposable, cotton	L	5 pcs
47512008*	SafeHandlingSheet, disposable, cotton	XL	5 pcs
47512009*	SafeHandlingSheet, disposable, cotton	XXL	5 pcs

* Without straps

SWL 300 kg / 660 lbs





TurningSling™

SystemRoMedic

For turning of users in bed

TurningSling is used, in combination with a two-point sling bar mounted on a stationary or mobile lift, for turning of users in bed. TurningSling also facilitates hygiene, treatment of sores and dressing.

- Ideal for coccyx, sacral, and spinal wound care and dressing changes
- Ideal for toileting patients who use a bedpan and for brief changes and patient hygiene including bed baths and linen changes
- Infection control (use with patients on contact, droplet and airborne precautions)
- Compatible with lifts and sling bars from other manufacturers

Item No.	Product	Size	Unit(s)
49592000	TurningSling, two loops, disposable, non-woven	one size	5 pcs
49592005	TurningSling, one loop, disposable, non-woven	one size	5 pcs

SWL 250 kg / 550 lbs



LimbSling™

SystemRoMedic

For easy lifting of an arm or leg

LimbSling is designed to comfortably and conveniently assist lifting of a user's leg or arm, for example, to gain better access for treatment and dressing sores. LimbSling is very practical for large and heavy users. LimbSling is applied by placing it under the limb to be lifted and connecting the strap loops to the sling bar.

- Ideal for wound care, limb holding, lymphedema, ROM therapy and turning
- Infection control (use with patients on contact, droplet and airborne precautions)
- Compatible with lifts and sling bars from other manufacturers

Item No.	Product	Size	Unit(s)
49292000	LimbSling, disposable, non-woven	one size	5 pcs

SWL 150 kg / 330 lbs



ThoraxSling™

SystemRoMedic

For sit-to-stand and standing training in combination with MiniLift200 and MiniLift125

ThoraxSling is a specially adapted sit-to-stand vest that is used in combination with Handicare's SystemRoMedic™ sit-to-stand lifts MiniLift200 and MiniLift125. It is easy for the caregiver to apply, while giving the user an added sense of security during rising and standing training. ThoraxSling also facilitates toileting.



Item No.	Product	Size	Unit(s)
45690004	ThoraxSling, with seat support, disposable, non-woven	S	5 pcs
45690006	ThoraxSling, with seat support, disposable, non-woven	M	5 pcs
45690007	ThoraxSling, with seat support, disposable, non-woven	L	5 pcs
45690008	ThoraxSling, with seat support, disposable, non-woven	XL	5 pcs

SWL 200 kg / 440 lbs



HighBackSling
SystemRoMedic™

BariSling
SystemRoMedic™

Standard model for users who need a lot of support

Handicare's disposable HighBackSling is designed for short-term use in most common lifting situations, e.g. sitting transfers between bed and wheelchair, and when the user needs extra head support. HighBackSling can also be used for lifting from a recumbent position on the floor.

Special model for extremely large and heavy users

Our disposable BariSling is a lifting sling that is designed for seated lifting of users weighing 200 kg/440lbs or more.

- Multiple uses during one user's hospital stay
- Sitting transfers (bed to chair, chair to chair, chair to bed, etc.)
- Infection control (use with patients on contact, droplet and airborne precautions)
- Compatible with lifts and sling bars from other manufacturers
- Bariatric option with up to 450 kg/1000 lbs capacity is available



Item No.	Product	Size	Unit(s)
40890004	HighBackSling, disposable, non-woven	S	5 pcs
40890006	HighBackSling, disposable, non-woven	M	5 pcs
40890007	HighBackSling, disposable, non-woven	L	5 pcs
40890008	HighBackSling, disposable, non-woven	XL	5 pcs
40890009	HighBackSling, disposable, non-woven	XXL	5 pcs

SWL 200 kg / 440 lbs, size XXL: SWL 300 kg / 660 lbs



47091010	BariSling, disposable, non-woven	XXXL	2 pcs
----------	----------------------------------	------	-------

SWL 450 kg / 1000 lbs



Disposable - but not for single use

All of our disposable transfer aids may be used for one user throughout their hospital stay. Disposable products facilitate infection control efforts. They are also a more cost effective option that eliminates the need for laundering as well as space required for large amounts of inventory.



HighBackSling/BariSling

StretcherSling

SafeHandlingSheet

TurningSling

ThoraxSling

LimbSling

EasyGlide

DisposableBelt

For updates and more information, please visit our website, www.handicare.com.



handicare

Maskinvägen 17,
972 54 Luleå, Sweden
Tel: +46 (0)8-557 62 200
Fax: +46 (0)8-557 62 299
www.handicare.com
info@handicare.se

SystemRoMedic™ is a complete solution for easy transfers. The concept encompasses four product areas that meet most needs in different transfer situations.

Transfer, between two locations. **Positioning**, within the same location.

Support, for mobility. **Lifting**, both manual and mechanical.